


# ENDOLINE® NASO-KIT

18

NASOBILIARY KIT

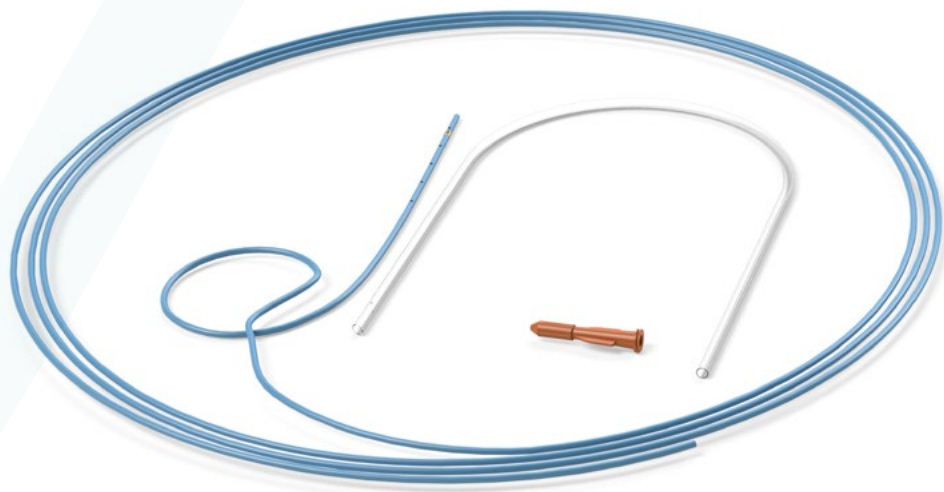

The ENDOLINE® NASO-KIT nasobiliary kit enables the aspiration of bile or pancreatic juices. This medical device is intended primarily for persons with obstructive jaundice, acute cholangitis or extrahepatic cholestasis.

## NASOBILIARY KIT COMPONENTS (components can be sold separately)

- A nasobiliary catheter for aspirating bile juice.
- A Luer-Lock connector mounted on the nasobiliary catheter enabling biliary or pancreatic juice aspiration. It connects the nasobiliary catheter to the full extension.
- A full extension to drain bile or pancreatic secretions in an external receiving bag or to rinse the catheter as needed.
- A nasal tube to retrieve the nasobiliary catheter in the larynx and to remove it through the nose.

## ADVANTAGES

- Axial and lateral perforations in the distal part allow for efficient drainage of the hepatobiliary tract.
- Easy to position catheter in the bile ducts thanks to the radiopaque distal end (with gold ring).

## AVAILABLE FORMATS (NASOBILIARY CATHETER)



CURVE 1



CURVE 2



CURVE 6

## ABOUT THIS MEDICAL DEVICE

1 <sup>st</sup> release date	1996
Single use	Yes
Class	Ila
Sterile	Yes
Sterilisation method	Ethylene oxide sterilisation
Contains latex	No
Lifetime	5 years
Certification	CE
Packaging	Box of 1 unit in individual sterile wrapping
Storage	In its original box in a cool and dry place

PRINCE MEDICAL SAS  
 ZA La Sente Du Moulin  
 64 rue Gutenberg - 60530 Ercuis - France  
[www.prince-medical.com](http://www.prince-medical.com)  
 Tel: +33 (0)3 44 26 54 76  
[prince-medical@omerin.com](mailto:prince-medical@omerin.com)

PRINCE MEDICAL  
 DEVICES

The information provided on this data sheet is indicative and subject to change without notice. PRINCE MEDICAL SAS shall not be held responsible for any incidents resulting from inappropriate use, failure to comply with the rules and regulations in force, or failure to follow the instructions for use provided. Medical devices sold by PRINCE MEDICAL SAS should only be used by a duly qualified practitioner or authorised person.

® Registered trademark of PRINCE MEDICAL SAS. Drawings and photos are not contractually binding. Reproduction is prohibited without the prior consent of PRINCE MEDICAL SAS.

## ENDOLINE® NASO-KIT

Reference	NASOBILIARY CATHETER		CONNECTOR LUER-LOCK	NASAL TUBE		EXTENSION
	Catheter diameter (Fr)	Catheter length (cm)	Connector Luer-Lock	Nasal tube Int. diam. / Ext. diam. (mm / mm)	Length of the nasal tube (cm)	Length of the extension (cm)
<b>CURVE 1</b>						
AS0540500	5	250	grey	2.0 / 3.0	50	135
AS0540700	7	250	orange	2.6 / 3.7	50	135
AS0541000	10	250	white	3.7 / 4.9	50	135
<b>CURVE 2</b>						
AS0560500	5	250	grey	2.0 / 3.0	50	135
AS0560700	7	250	orange	2.6 / 3.7	50	135
<b>CURVE 6</b>						
AS0580500	5	250	grey	2.0 / 3.0	50	135
AS0580700	7	250	orange	2.6 / 3.7	50	135

PRINCE MEDICAL SAS  
 ZA La Sente Du Moulin  
 64 rue Gutenberg - 60530 Ercuis - France  
[www.prince-medical.com](http://www.prince-medical.com)  
 Tel: +33 (0)3 44 26 54 76  
[prince-medical@omerin.com](mailto:prince-medical@omerin.com)

The information provided on this data sheet is indicative and subject to change without notice. PRINCE MEDICAL SAS shall not be held responsible for any incidents resulting from inappropriate use, failure to comply with the rules and regulations in force, or failure to follow the instructions for use provided. Medical devices sold by PRINCE MEDICAL SAS should only be used by a duly qualified practitioner or authorised person.

® Registered trademark of PRINCE MEDICAL SAS. Drawings and photos are not contractually binding. Reproduction is prohibited without the prior consent of PRINCE MEDICAL SAS.